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USER MANUAL

BeautyDock+™ | ALL-IN-ONE SKIN DEVICE

STRAFFEND
FIRMING

REGENERIEREND
REVITALIZING

STÄRKEND
STRENGTHENING

ABSCHWELLEND
DEPUFFING

DEFINIEREND
CONTOURING

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
Enjoy every moment with your
BeautyDock+™ | All-in-One Skin Device with
EMS and Red Light from Rosental Organics!

Before using your new tool, please read these instructions carefully and keep them for future reference.

To ensure that you use this product correctly, please read all safety instructions carefully before using it. The safety instructions will help you to use this product correctly to avoid personal injury or property damage to yourself or others. You will find the safety instructions below. Please observe the relevant regulations!


PRODUCT OVERVIEW

INCLUDED IN THE BOX

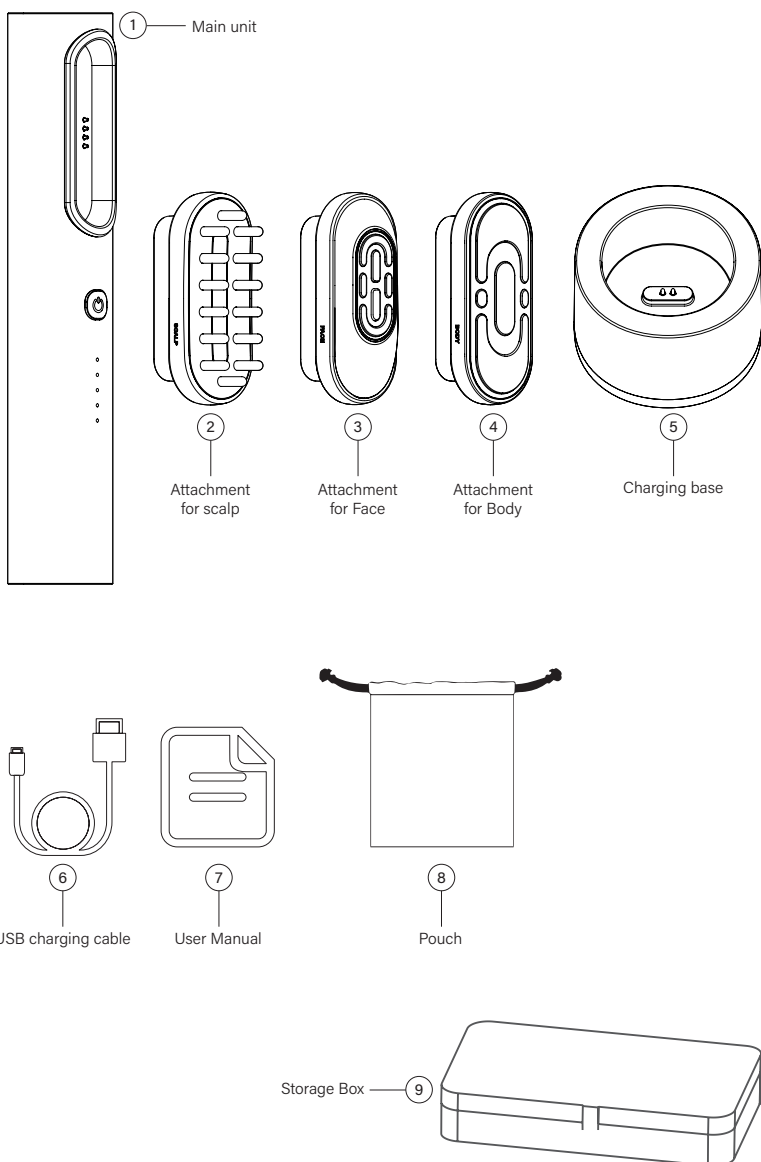PRODUCT OVERVIEW

The LED indicator displays the battery level when the battery is low or while the device is charging.

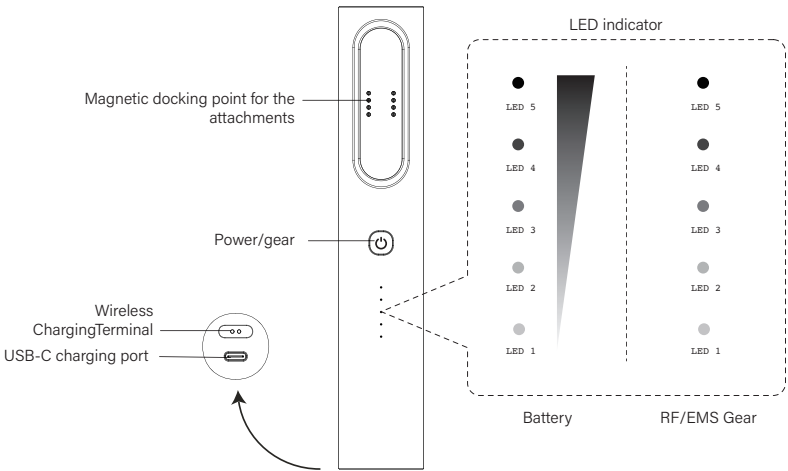
PRODUCT OVERVIEW

ATTACHMENTS

The function mode switches automatically when you change the attachment.


SCALP ATTACHMENT:

Ideal for relaxing and stimulating scalp care.

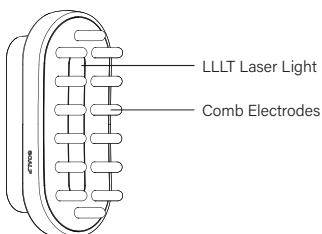
FACE ATTACHMENT:

Helps to firm the skin and define facial contours.


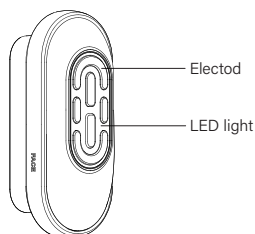
BODY ATTACHMENT:

For targeted care and firming of body areas.

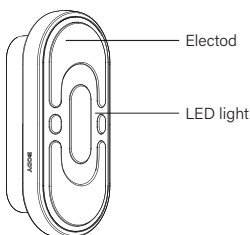
Carefully insert the attachment into the corresponding slots in the device until it clicks securely into place. When inserting the attachments, make sure the label „FACE/BODY/SCALP“ is facing downward. To replace an attachment, simply remove the current one and attach the desired one. The device will automatically recognize the new attachment and switch to the corresponding mode.

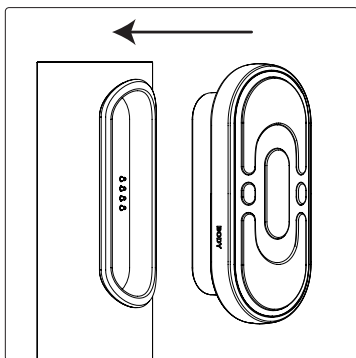
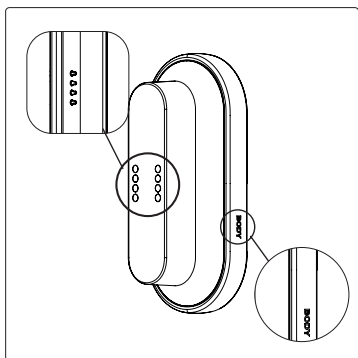
PRODUCT OVERVIEW

CHARGING BASE

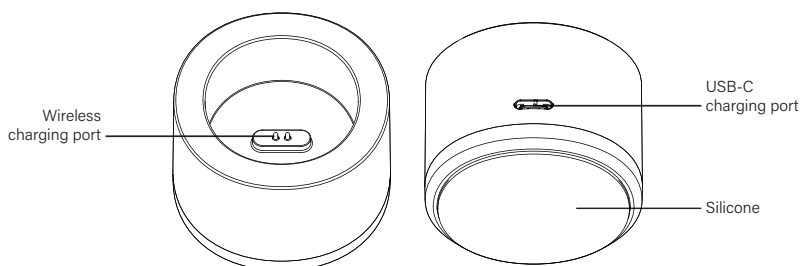The charging base is not waterproof. Please keep it in a dry place. Never use or charge the device in wet conditions or in the bathroom.

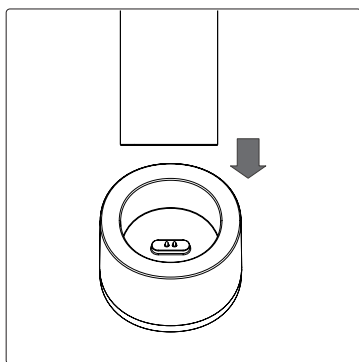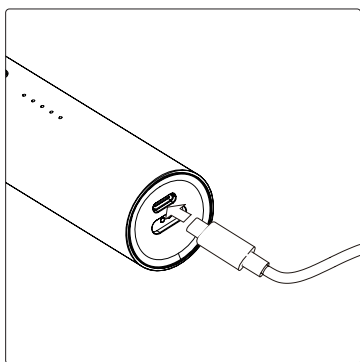CONTACT CHARGING:

Fully insert the USB-C charging cable connector into the charging port of the base, then place the main unit directly into the charging base. Ensure the device is completely seated in the base. The LED indicator light will turn on, indicating the device is successfully powered.

PLUG-IN CHARGING:

When the USB-C charging cable connector is fully inserted into the charging port of the main unit, the LED indicator light will turn on to indicate the device is powered on and charging.

When all 5 LED indicator lights remain steadily lit, the battery is fully charged.

When the battery is low, the device will vibrate 5 times and the 5 LED lights will flash 5 times.

The device cannot be used while charging.

The product should be charged lying flat on a surface or stored in the charging base.

MODES

The device automatically detects the attached attachment. As soon as you attach your desired attachment (Face, Body, or Scalp), the corresponding mode is automatically activated.

MODES	APPLICATION AREA	USAGE TIME	USAGE RECOM-MENDATIONS	INTENSITY LEVELS
FACE	Face and neck	10 Minutes	3-4 times per week for no more than 20 minutes each time	5 Intensity levels
BODY	Body (arms, thighs, back, abdomen, etc.)	10 Minutes	3-4 times per week for no more than 20 minutes each time	
SCALP	Scalp	10 Minutes	Use once a day for no more than 30 minutes	

NOTICE: BEFORE USE

- The sensation of warmth, mild electrical pulses or vibrations during the use of the **BeautyDock+™** is normal and intentional. These technologies activate the skin cells, promote collagen production, and support scalp health.
- Before each use, check the **BeautyDock+™**. The device should not be used if it shows obvious signs of damage. If the device is damaged due to shocks or falls, stop using it immediately. If the tool is damaged due to shocks or falls, please do not continue to use it.
- The device cannot be operated during the charging process. Only use the included charging base and cable to charge the **BeautyDock+™**. The controller's indicator light blinks during charging. Do not twist, bend, forcibly pull out, wrap, or otherwise damage the charging cable, as this could lead to fire or short circuit. If the charging cable is damaged or becomes hot, stop using it immediately.
- Do not store the device near fire sources or at very high ambient temperatures. Keep the tool away from strong electromagnetic fields, such as induction fields, during use. Avoid contact with strong electromagnetic fields, such as those produced by induction devices.

NOTICE: BEFORE USE

- Direct applications to the eyes, eyeballs, wounds, plastic parts, neck, larynx, or over the heart should be avoided. To protect your eyes: Do not look directly into the LED light. Ensure that the LED lights do not shine into your eyes during use. If improper use or use of the product leads to unwanted skin reactions or other discomfort, it is recommended to seek medical advice.
- The LED red light therapy works in the dermis, where the energy is precisely confined to a safe area and does not harm the skin.
- The device should not be immersed in water or other liquids.
- This product is not suitable for children, pregnant or breastfeeding women, individuals with allergies, skin injuries, eczema, swelling, skin diseases, atopic dermatitis, or sensitive skin. Likewise, people with mental illnesses, pacemakers, or other devices, individuals with fever, heart diseases, severe or acute illnesses, malignant tumors, tuberculosis, or blood disorders should not use this device.
- The device is not a toy. Keep it out of the reach of children (0-8 years old) and ensure that they do not perform any maintenance or cleaning on the device.
- Do not alter, disassemble, or repair the device yourself. For repairs, please contact the distributor or manufacturer directly. Do not use accessories that are not part of the supplied contents.
- This product is an electronic device for home use, not a medical device, and is not intended for diagnostic or therapeutic purposes.


APPLICATION: BEFORE FIRST USE

Perform a sensitivity test before your first application:

1. Place the **BeautyDock+™** on your inner arm.
2. Press and hold the ON/OFF button until the lights turn on, and allow the device to complete the full treatment program.
3. Make sure to test the device at the lowest setting.
4. If your skin feels too hot or uncomfortable at any point, turn off the device immediately and remove it.
5. Wait at least 6 hours after the treatment. If redness or itching occurs on your inner arm, the **BeautyDock+™** may not be suitable for you. Please consult your doctor before use.


This product conducts electricity by making contact with the skin using both positive and negative electrodes simultaneously. Effective stimulation occurs when the electrode needle or the entire electrode touches the skin.

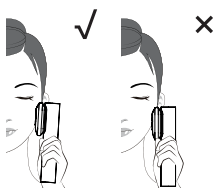
When using the Face Attachment:

Try to apply the treatment surface as fully as possible to the skin, and avoid floating it during treatment.

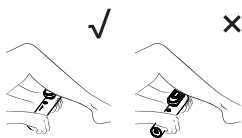When using the Body Attachment:

Try to keep the treatment surface as close to the skin as possible, and avoid floating.

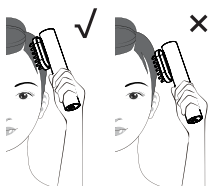
When using the Scalp Attachment:

Try to fit the entire side of the brush head to the scalp and avoid floating during treatment.


MODE 1:

FACE

1. After attaching the face attachment, press and hold the power button for 1.5 seconds to turn on the device. The red and infrared lights will turn on.
2. Shortly press the power button to start the treatment. It is recommended to start with intensity 1 and adjust the level according to your skin's sensitivity. The intensity cycles in the order: „1-2-3-4-5-1-2...“.

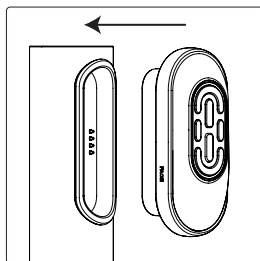3. Facial Massage Techniques: Follow the illustrated steps for massage. Steps ① to ④ form one complete cycle. Repeat the steps 2 to 3 times on each side of the face. The total treatment time for both sides of the face is about 10 minutes. When the treatment ends, the device will vibrate twice and automatically turn off.

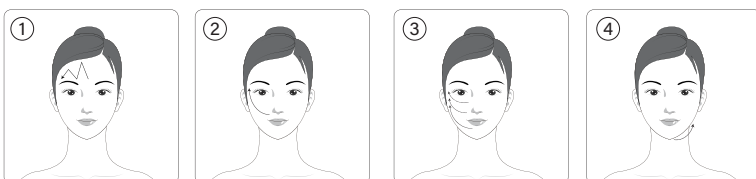
- ① **Forehead:**
Move in a Z-shaped motion from between the eyebrows toward the temples.
- ② **Cheeks:**
Glide from the corners of the mouth toward the ears.
- ③ **Smile Lines:**
Glide upwards along the smile lines, then move toward the temples.
- ④ **Jawline:**
Move from the chin toward the earlobes.




4. To power off, press and hold the power button for 1.5 seconds again, then remove the face attachment.

MODE 2: BODY

1. After attaching the body attachment, press and hold the power button for 1.5 seconds to turn on the device. It will vibrate briefly, and the red LED light will turn on to indicate startup.
2. Short-press the power button to start the treatment. It's recommended to begin with intensity 1 and adjust the intensity according to your skin's sensitivity. The intensity levels cycle in the following order: "1-2-3-4-5-1-2..."

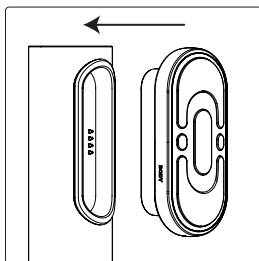
3. Point-Press Massage Technique: Using the shoulders and neck as an example (see illustration), follow these steps: Steps ① to ③ form one complete massage cycle. Hold each pressure point for about 8 seconds, and repeat the full cycle 2–3 times per side. The total recommended treatment time for both sides is around 10 minutes. After the session, the device will vibrate twice and automatically power off.

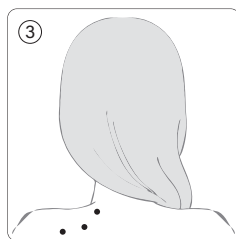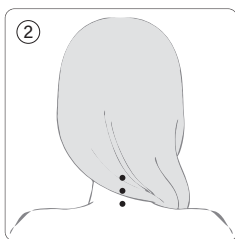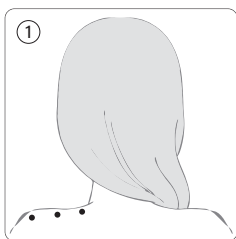- ① **Shoulders:**
Start from the base of the neck and press outward toward the arms.
- ② **Neck:**
Apply steady pressure downward along the nape of the neck.
- ③ **Upper back (trapezius area):**
Start at the neck base and press outward toward the trapezius muscle.

MODUS 2: KÖRPER


4. Sliding Massage Technique: For other body areas, you can use the following techniques (see illustration):
-



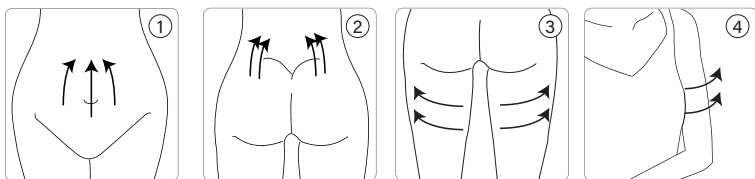
- ① **Abdomen:**
Glide slowly from bottom to top to firm the area.
- ② **Glutes:**
Glide slowly upward to ease nerve tension.
- ③ **Legs:**
Glide slowly from the inside outward and from bottom to top to tone the legs and relieve muscle tension.
- ④ **Arms:**
Glide slowly from the inside outward and from bottom to top to shape the arms and ease muscle tension.
-
5. After completing your massage, press and hold the power button for 1.5 seconds to power off, then remove the body attachment.

MODE 3: SCALP

1. After attaching the scalp attachment, press and hold the power button for 1.5 seconds to turn on the device. It will vibrate briefly, and the red LED light will turn on to indicate startup.
2. Short-press the power button to start the treatment. It's recommended to begin with intensity 1 and adjust the intensity according to your skin's sensitivity. The intensity levels cycle in the following order: "1-2-3-4-5-1-2..."

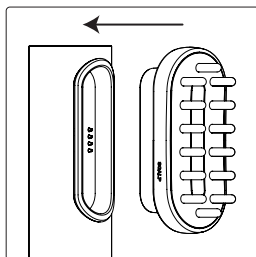3. Scalp Massage Technique: As shown in the illustration, follow these three steps for a full massage cycle: Each section should take about 8–10 seconds. Repeat the entire cycle 2–3 times per side. The total care time for the whole scalp is about 10 minutes. Once finished, the device will vibrate twice and automatically turn off.

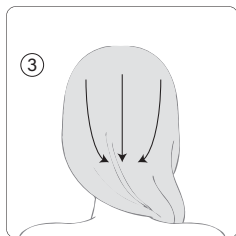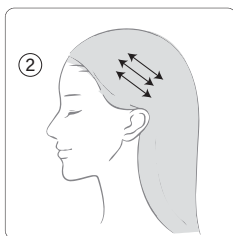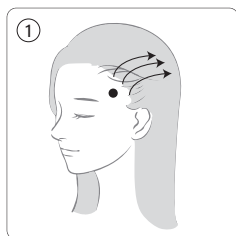
- ① Slowly comb backward from the temples.
- ② Move slowly in a gentle zigzag pattern from the temples to the area above the ears.
- ③ Glide slowly from the crown of your head down to the back of the head, keeping the brush close to the scalp.

4. After the massage, press and hold the power button for 1.5 seconds to turn off the device and remove the scalp attachment.

The device will automatically power off if left idle for 3 minutes in standby mode.

Tips:

1. Do not look directly into the red light for extended periods to avoid eye damage.
2. Do not use the metal comb teeth near the eye area to prevent eye injury.
3. The recommended treatment duration per area is at least 10 minutes, up to a maximum of 30 minutes.
4. After use, gently wipe the device with a damp paper towel or soft cloth.

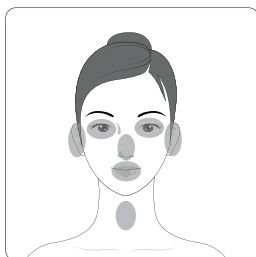CLEANING AND STORAGE

To ensure optimal performance and extend the lifespan of the device, always place it back into the charging base after each use. This helps prevent dust buildup and accidental damage.

To clean the device, gently wipe the main unit with a soft, damp cloth.

The attachments can be rinsed with water if necessary, then dried with a soft cloth.

Please note: The base is not waterproof. Before cleaning the base, unplug the power supply and wipe it gently with a soft, damp cloth.

BATTERY

The **BeautyDock+™** is powered by an integrated rechargeable lithium polymer battery, which must be properly recycled at the end of its life cycle. Please always follow the correct instructions below when charging the battery:

1. Fully connect the supplied Type-C charging cable to the charging port of the device. During the charging process, the intensity indicator will start flashing, starting with the first level. The battery is fully charged as soon as all three intensity/battery status indicators are permanently lit. A full charge takes about 2 hours.
2. Charge the device in good time if the battery is low. This is indicated by all mode and intensity indicators flashing five times and five signal tones.
3. Disconnect the device from the power source after fully charging to avoid overcharging, which can affect battery life.
4. The device cannot be operated during the charging process.
5. Only charge the device at ambient temperatures between 0°C/32°F and 40°C/104°F. An ambient temperature of 10°C-30°C is recommended. Note that the charging process may be slower in an environment below 10°C or when charging via a computer.
6. Lay the device on its side while charging and avoid charging in humid environments.
7. If it is not used for a longer period of time, it is recommended to fully charge it every 3 months to prevent damage to the battery.

CAUTION!

1. The following persons should not use this product:
 - People under 18 years old
 - Pregnant or breastfeeding women
 - People with allergies
 - People with skin injuries, eczema, swelling, skin diseases, atopic dermatitis or sensitive skin
 - People with any pre-existing inflammation
 - People with mental illnesses and impaired mental abilities
 - People with a pacemaker or other devices
 - People with fever and other physical abnormalities, heart disease, serious illnesses, acute illnesses, malignant tumours, tuberculosis, blood and blood pressure disorders and other symptoms
 - People with epilepsy or a history of skin cancer (consultation with a doctor is required before use)
 - People with any existing inflammation or pre-existing conditions
2. Do not use it on eyes, eyeballs, wounds, plastic parts, anterior neck or larynx, tumours, eczema, over the heart, etc.
3. Do not use the tool in case of any pre-existing inflammations or other illnesses.
4. Avoid using the device after recent injuries, surgeries, or aesthetic treatments, including neurotoxin treatments, dermal fillers, microneedling, laser treatments, and chemical peels, until the skin has fully healed.
5. After aesthetic treatments such as fillers, hyaluronic acid, Botox, etc., consult your doctor before use.
6. If you experience discomfort during use, stop immediately and seek medical advice.
7. Please do not exceed the recommended usage time.
8. The different modes of the **BeautyDock+**™ can generally be used in succession, but ensure that the total application time is not exceeded.
9. If you have sensitive skin, consult a dermatologist before the first use.
10. Individuals who are allergic to light or have light allergies due to medication use should not use this device
11. Do not look directly into the LED light for an extended period, as this may cause eye damage.
12. If redness or swelling of the skin lasts longer than 2 hours, discontinue use and contact a medical professional.
13. Do not use the device on open wounds.

CAUTION!

14. Do not use the device on painful, raised, crusty, or bleeding skin. This could be a sign of a more serious condition and should be checked by your doctor.
15. Use the **BeautyDock+™** only for the conditions specified in the application guidelines. The device has been tested exclusively for these listed conditions, and the risks for other conditions are unknown.
16. Your skin may become more sensitive to sunlight for two or three days—wear a hat or apply sunscreen.
17. Avoid contact with water during use.
18. Make sure your hands are dry when connecting or disconnecting the charging station.
19. Do not store the tool near a fire source, heat it, or use it in very high ambient temperatures. Store it in a normally tempered, dry, and well-ventilated place.
20. If the tool is damaged due to shocks or falls, do not continue using it.
21. Do not modify, disassemble, or repair the tool yourself. For repairs, please contact the distributor or manufacturer directly.
22. Prevent the internal battery from coming into contact with skin, eyes, or mouth. If exposed to battery fluids, wear protective gloves and seek medical assistance in case of skin or eye contact.
23. Do not use accessories that are not included in the original package.
24. Do not twist, bend, forcibly pull out, wrap, or otherwise damage the charging cable, as this could cause fire or a short circuit. If the charging cable is damaged or becomes hot, stop using it immediately.
25. Avoid operating the device under blankets or pillows to prevent excessive heat buildup and potential fire hazards.

CARE AND MAINTENANCE

1. Fully charge the device before using it for the first time.
2. Avoid washing the device with a strong jet of water or immersing it directly in water.
3. Do not use steel wool, alkaline cleaning agents, abrasive cleaning agents, alcohol, thinner, petrol, nail polish remover or other corrosive liquids for cleaning.
4. The device uses lithium batteries. Please do not place or store it in a humid or high-temperature environment for long periods of time.
5. If the tool is not used for a long time, it should be switched on regularly and recharged every three months to check its functionality. Before prolonged periods of non-use, clean the tool and store it in a closed and dry place.
6. Do not open the tool, batteries or charger as this may result in a risk of electric shock, fire or exposure to battery chemicals. Disassembly or incorrect reassembly will invalidate the warranty.
7. Only use the specified cleaning products to avoid damage.

DISPOSAL

Before disposal, dismantle the tool properly and remove the battery from the remote control. First, turn off the power supply.

Carefully disconnect the battery connector from the circuit board and remove the battery.

Dispose of the battery separately from the housing properly at a recycling station according to local guidelines.

This symbol indicates that batteries and accumulators must not be disposed of with regular household waste. You are legally obligated to return them free of charge to a publicly authorized collection point or to collection points set up by distributors for recycling at the end of their lifespan.

Batteries and accumulators that are not firmly enclosed in the device and can be removed without destruction must be separated and disposed of in a discharged state according to the designated disposal method. The recycling and proper disposal of old devices and batteries are important contributions to protecting our environment. Improper disposal can lead to toxic substances entering the environment, which can have harmful effects on humans, animals, and plants.

Below are some best practices and recommendations to extend the lifespan of your batteries and accumulators and enable their reuse:

- Charge your batteries properly and fully to maximize their lifespan. Fully discharge them with a suitable charger if necessary before recharging.
- Consider using rechargeable batteries. These can be reused and recharged multiple times, reducing the amount of battery waste.
- Always use the correct type of batteries for your devices. Improper usage can shorten battery life and may have harmful effects.

Always follow safety instructions when handling used batteries. Lithium-containing used batteries (Li = Lithium) pose a high fire hazard. Therefore, special attention must be paid to the proper disposal of lithium-containing used batteries and accumulators. Improper disposal can also lead to internal and external short circuits due to thermal effects (heat) or mechanical damage. A short circuit can result in a fire or explosion and have serious consequences for humans and the environment. To prevent a short circuit, tape over the terminals or contacts of the battery before disposal.

TROUBLESHOOTING

If a malfunction occurs during operation, please contact the customer support directly via **contact@rosental.de** and do not disassemble and attempt to repair it by yourself.

WARRANTY

The following conditions, which describe the requirements and scope of our warranty service for the product **BeautyDock+™** | All-in-One Skin Device with EMS and Red Light by Rosental Organics, do not restrict your statutory rights in the event of product defects. You may exercise these statutory rights free of charge, regardless of this warranty. In addition, we grant end users of the product a warranty under the following conditions:

1. We will, free of charge and in accordance with the following conditions (clauses 2 to 9), remedy defects in the product that can be proven to be due to material and/or manufacturing faults, provided such defects are reported immediately upon discovery and within 12 months of delivery to the original end customer.
2. Product defects or damage caused by human error, misuse, or improper use (contrary to this instruction manual) are excluded from the warranty.
3. Severe impacts or collisions, excessively high or low power supply voltages, as well as product defects or damage caused by unavoidable circumstances, are excluded from the warranty.
4. The warranty service will be provided at our discretion either by (i) repairing the defective parts of the product free of charge or replacing them with flawless parts (repair), or (ii) replacing the product with a defect-free one (replacement). Replaced parts or devices become our property.
5. You can claim rights under this warranty by contacting customer service directly at [**contact@rosental.de**](mailto:contact@rosental.de) (with presentation of the original invoice showing the purchase date).
6. For the above-mentioned cases not covered by the warranty, we will decide whether a repair is possible, offer you the repair, and then invoice you accordingly for parts and repair costs.
7. This warranty applies only to products purchased and used by consumers for private purposes, not for industrial or commercial use.
8. This warranty applies to devices purchased and operated within the European Union.
9. Warranty services do not extend the original warranty period nor initiate a new period. The warranty period for replaced or repaired parts ends at the same time as the warranty period for the entire product.
10. This warranty is governed by the laws of the Federal Republic of Germany.
11. Warranty provider (Manufacturer & Distributor):
Rosental Organics GmbH
Rosa-Luxemburg-Str. 14,
10178 Berlin, Germany

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