

# SONIC KIDS TOOTHBRUSH

Set

 **geokids**  
by georganics



## User Guide




- 1.** Straight from the box, start using your toothbrush until the battery is completely discharged of power. Then fully charge your toothbrush for 12 hours




Please read the safety introduction carefully before use

- 2.** Align the front of the brush head with the front of the toothbrush handle. Insert the brush head into the metal shaft.

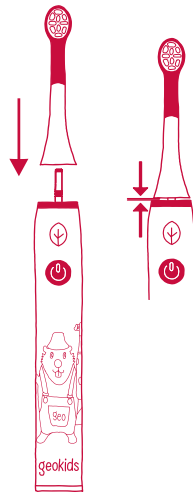

Please read the safety introduction carefully before use

- 3.** Apply toothpaste onto bristles, horizontally.

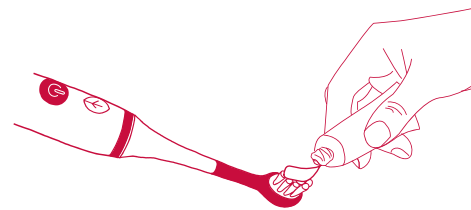


Please read the safety introduction carefully before use

- 4.** Position the bristles between teeth and gum line gently, on 45 degree angle, applying low to moderate pressure.


### WARNING

Children must be supervised to ensure proper usage. The product is not a toy.

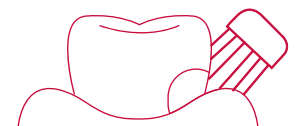


Please read the safety introduction carefully before use

- 5.** Turn the toothbrush ON and clean both lower and upper arch spending at least 30 seconds on each half section.

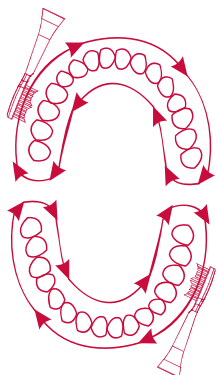


Please read the safety introduction carefully before use

- 6.** Turn the toothbrush OFF while still in your mouth.

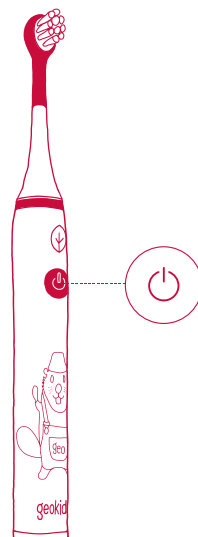

Please read the safety introduction carefully before use


- 7.** Using water, rinse the brush head after each use and wipe away water from the handle by using a damp cloth.

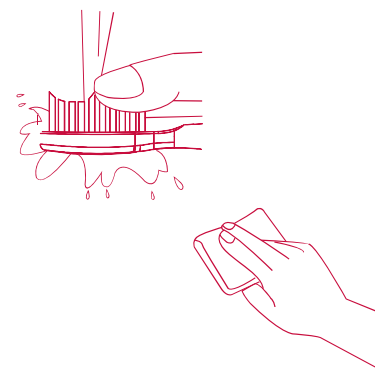

Please read the safety introduction carefully before use

- 8.** Repeat twice a day as recommended by dentist.

### BRUSH TWICE A DAY

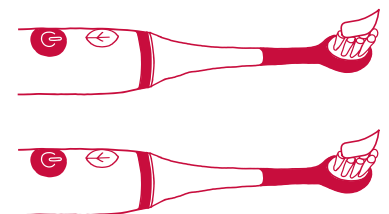

Please read the safety introduction carefully before use